

# Nio

## User Guide



Supported displays:

E-3620, MDNG-5121, MDNG-2121,  
MDNC-3121, MDNG-6121, MDNC-2121

**BARCO**







Visibly yours

(This page intentionally left blank.)

# Welcome!

Nio is Barco's flexible, industry-standard display solution for diagnostic imaging. Presenting an exceptional combination of innovation and performance, Barco's next-generation Nio display system improves your diagnostic reading routine for a multitude of medical imaging applications.

Following symbols may be used throughout this guide:

	Warning: Risk of injury to human beings		Caution: Risk of damage to the product
	Important notice or remark		Note
	Hint, tip		Additional information



Read all the safety information starting on page 43 before operating your BARCO Nio Display.

(This page intentionally left blank.)

# Table of contents

<b>What's in the box.....</b>	<b>7</b>
<b>Parts, controls and connectors .....</b>	<b>8</b>
Front view .....	8
Rear view .....	9
<b>Display installation .....</b>	<b>12</b>
Unlocking the height mechanism.....	12
Adjusting the display position .....	13
Connecting the signal cables.....	16
Routing the cables.....	17
Re-attaching the covers .....	17
VESA-mount installation.....	18
<b>Daily operation .....</b>	<b>20</b>
Recommendations for daily operation .....	20
Stand-by switching .....	22
Locking and unlocking user controls.....	23
Bringing up the OSD menu.....	25
Navigating through the OSD menus.....	26
Making selections in the OSD menus.....	26
Changing values in the OSD menus .....	27
<b>Advanced operation .....</b>	<b>28</b>
On-screen display (OSD).....	28
Input Selection menu* .....	28
Luminance and color menu .....	29
Display Function menu* .....	30
Settings menu .....	33
Information menu.....	34
Changing Display Functions.....	35

Concepts.....	35
Display function selection .....	36
ALC & DICOM Options.....	37
Switching DPMS on/off.....	40
<b>Cleaning your display .....</b>	<b>41</b>
Front glass .....	41
Cabinet.....	42
<b>Important information .....</b>	<b>43</b>
Safety information .....	43
Environmental information.....	46
Regulatory compliance information.....	49
Explanation of symbols .....	50
Legal disclaimer .....	52
Technical specifications .....	53
Warranty statement.....	71

# What's in the box

Your BARCO NIO Display comes with:

- this BARCO NIO Display guide
- a system CD
- a DVI cable
- a DisplayPort cable (MDNC-2121 only)
- a USB cable
- a set of AC power cords
- an external power supply

If you ordered a Barco display controller, it's also in the box together with its accessories. A dedicated user guide is available on the system CD.



Keep your original packaging. It is designed for this display and is the ideal protection during transport.

# Parts, controls and connectors

## Front view

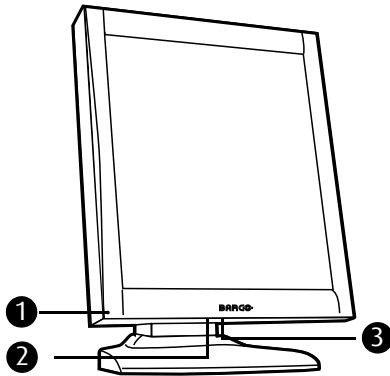


Figure 1: Front side

---

1 Power LED

2 USB downstream port

---

3 Control wheel

---



# Rear view

## REAR VIEW OF E-3620, MDNG-5121

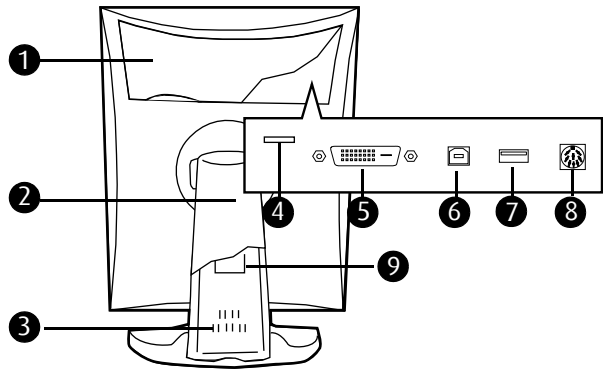


Figure 2: Rear side

1	Connector compartment cover	2	Tilt & swivel foot cover
3	Tilt & swivel foot	4	Slot for security cable
5	DVI (digital) video input	6	USB upstream port
7	USB downstream port	8	DC power input
9	Tilt & swivel foot clip		

# REAR VIEW OF MDNC-3121, MDNG-2121 AND MDNG-6121

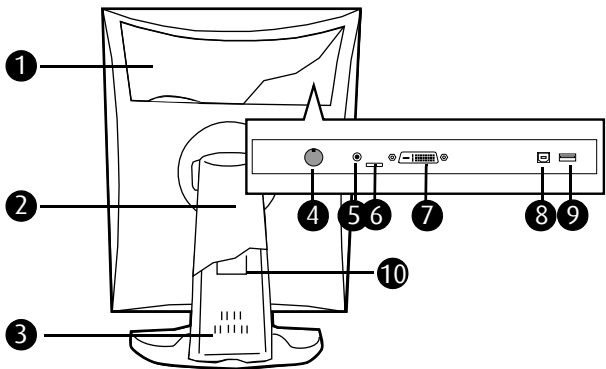


Figure 3: Rear side

1	Connector compartment cover	2	Tilt & swivel foot cover
3	Tilt & swivel foot	4	DC power input
5	+5 Vdc	6	Slot for security cable
7	DVI (digital) video input	8	USB upstream port
9	USB downstream port	10	Tilt & swivel foot clip

# REAR VIEW OF MDNC-2121

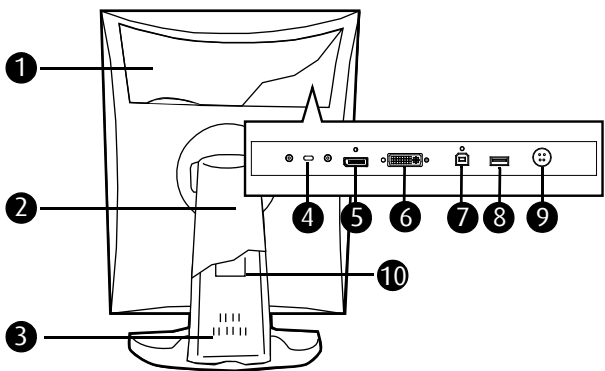


Figure 4: Rear side

1	Connector compartment cover	2	Tilt & swivel foot cover
3	Tilt & swivel foot	4	Slot for security cable
5	DisplayPort video input	6	DVI (digital) video input
7	USB upstream port	8	USB downstream port
9	DC power input	10	Tilt & swivel foot clip

# Display installation



Prior to installing your BARCO Nio Display and connecting all necessary cables, make sure to have a suitable display controller physically installed in your computer. If you are using a Barco display controller, please consult the dedicated user guide available on the system CD.

## Unlocking the height mechanism

In the factory, the height-positioning system in the display foot is blocked with a red clip to prevent damage during transportation. Before installing the display, you must remove this clip.

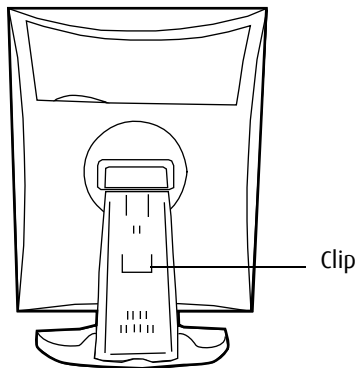


Figure 5: Location of the clip

### TO REMOVE THE CLIP:

- 1 Position the display with its rear side facing you.
- 2 Pull the red clip out of the fixation holes in the foot.
- 3 Keep the clip in case the display needs to be shipped later.

## Adjusting the display position

You can change the orientation of the panel at any time, but it is more convenient to select landscape or portrait orientation before connecting the cables..

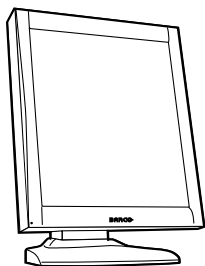


Figure 6: Portrait orientation

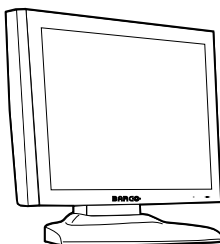


Figure 7: Landscape orientation

## TO CHANGE THE PANEL ORIENTATION

- 1 Stand at the front side of the panel and take the panel at both sides.
- 2 Very important: Til the panel before changing the orientation.



Should you change the panel orientation without tilting it first, you might irreversibly damage the tilt & swivel mechanism.

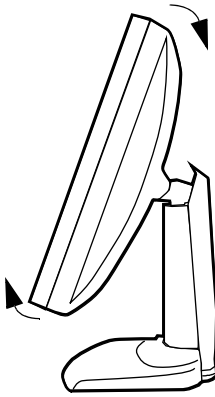


Figure 8: Tilt the panel before rotating

- 3 To change from portrait to landscape, turn the panel counterclockwise.

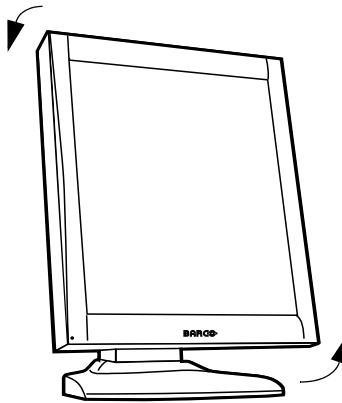


Figure 9: To rotate the panel from portrait to landscape

- 4 To change from landscape to portrait, turn clockwise..



If, after installing the display or the system, you change the panel orientation while an image is on the screen, the result depends on the graphic board and the resolution of the image. In some cases the image will be rotated automatically, in other cases it will not be rotated (e.g., when pixels would be lost after rotation). If necessary, change the image resolution in the display control panel and restart the system after changing the orientation.

# Connecting the signal cables

## TO CONNECT THE SIGNAL CABLES TO THE DISPLAY:

To get access to the connectors, remove the connector compartment cover by pulling down the 2 clips at the top of the cover.

The location of the connectors depends on the display type:

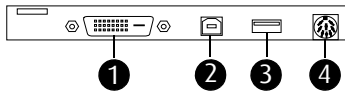


Figure 10: Inputs E-3620, MDNG-5121

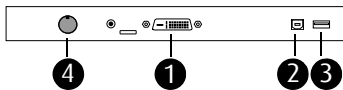


Figure 11: Inputs MDNC-3121, MDNG-2121 and MDNG-6121

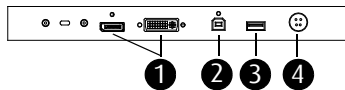


Figure 12: Inputs MDNC-2121

- Connect one end of the DVI/DisplayPort cable to the DVI/DisplayPort input of the display (1). Connect the other end of the DVI/DisplayPort cable to the DVI/DisplayPort connector of the display controller board.
- Connect a PC USB downstream connector to the display's USB upstream connector by means of a USB cable (2).
- Connect any USB device to one of the display's USB downstream connectors (3)
- Connect the DC power input (4) of the display to the external DC power supply. Connect the other end of the external DC power supply to a **grounded** power outlet by means of the proper power cord delivered in the packaging.



## Routing the cables

- Bind the cables in the connector compartment together with the cable tie inside the connector compartment.
- Put the connector compartment cover back on the display. Pay attention that the signal cables are positioned under the bulge in the cover.
- Push the cables into the clips on the rear of the tilt & swivel foot.
- Bind the cables together above and under the foot, by means of the 2 velcro strips attached to the inside of the foot cover (packed inside the accessory box).
- At last, put the foot cover back in place.

## Re-attaching the covers

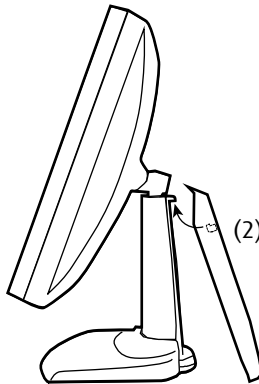


Figure 13: Installing the cover

- 1** The foot cover is packaged separately in the accessory box. Unpack the foot cover.
- 2** Push the upper side of the cover onto the foot, so that the hooks inside the cover are positioned right under the bulges at the rear of the foot.

- 3** Slide the cover upward while moving the lower side of the cover towards the foot.
- 4** Press the cover to the foot so that it makes a clicking sound.

## VESA-mount installation

The panel, standard attached to the tilt & swivel foot, is compatible with the VESA 100 mm standard. So it can be used with an arm stand according to the VESA 100 mm standard.

Therefore, the tilt & swivel foot must be removed from the panel..



- Use an arm that is approved by VESA (according to the VESA 100 mm standard).
- Use an arm that can support a weight of at least 13 kg (28.66 lbs).

### TO ATTACH THE DISPLAY TO AN ARM STAND:

- 1** Put the display face down on a clean and soft surface. Be careful not to damage the panel screen.
- 2** Remove the tilt & swivel foot cover.

- 3 Remove the small screw (A) fixing the small plastic cover on top of the foot. Next, remove the small cover itself.

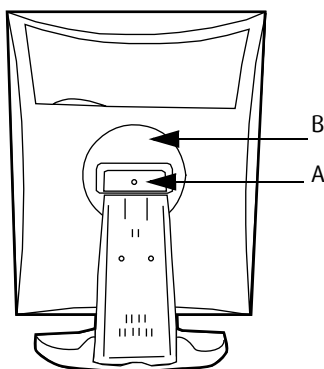


Figure 14: Display with tilt & swivel foot cover removed

- 4 Unscrew the 2 screws fixing the round plastic cover (B).
- 5 Lift up the round plastic cover.
- 6 Remove the four screws fixing the foot while supporting the foot.
- 7 Attach the arm stand **firmly** to the panel using 4 screws M4 x 8 mm.

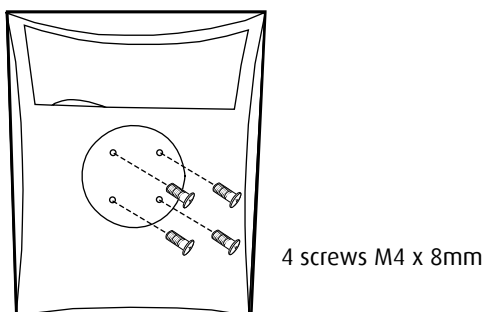


Figure 15: Location of the screws

# Daily operation

## Recommendations for daily operation

### OPTIMIZE THE LIFETIME OF YOUR DISPLAY

Enabling the Display Power Management System (DPMS) of your display will optimize its diagnostic lifetime by automatically switching off the backlight when the display is not used for a specified period of time. By default, DPMS is enabled on your display, but it also needs to be activated on your workstation. To do this, go to “Power Options Properties” in the “Control Panel”.



Barco recommends setting DPMS activation after 20 minutes of non-usage.

### USE A SCREEN SAVER TO AVOID IMAGE RETENTION

Prolonged operation of an LCD with the same content on the same screen area may result in a form of image retention.

You can avoid or significantly reduce the occurrence of this phenomenon by using a screen saver. You can activate a screen saver in the “Display properties” window of your workstation.



Barco recommends setting screen saver activation after 5 minutes of non-usage. A good screen saver displays moving content.

In case you are working with the same image or an application with static image elements for several hours continuously (so that the screen saver is not activated), change the image content regularly to avoid image retention of the static elements.

### UNDERSTAND PIXEL TECHNOLOGY

LCD displays use technology based on pixels. As a normal tolerance in the manufacturing of the LCD, a limited number of these pixels may remain either dark or permanently lit, without affecting the diagnostic

performance of the product. To ensure optimal product quality, Barco applies strict selection criteria for its LCD panels.



To learn more about LCD technology and missing pixels, consult the dedicated white papers available at [www.barco.com/medical](http://www.barco.com/medical).

## ENHANCE USER COMFORT

Every Barco multi-head display system is color matched with the highest specification in the market.



Barco recommends keeping color-matched displays together. Furthermore, it is important to use all displays of a multi-head configuration at the same rate to preserve color matching throughout the economic lifetime of the system.

## MAXIMIZE QUALITY ASSURANCE

The 'MediCal QAWeb' system offers online service for high-grade Quality Assurance, providing maximum diagnostic confidence and uptime.



Barco recommends to install MediCal QAWeb Agent and apply the default QAWeb policy at least. This policy includes calibration on regular intervals. Upgrading to MediCal QAWeb Server offers even more possibilities. Learn more and sign up for the free MediCal QAWeb Essential level at [www.barco.com/medical/qa](http://www.barco.com/medical/qa)

# Stand-by switching

## ON E-3620 AND MDNG-5121 DISPLAYS:

When the display is on and no on-screen display is visible, push and hold the control wheel at the front for a few seconds to switch the display in stand-by. The LED turns orange.

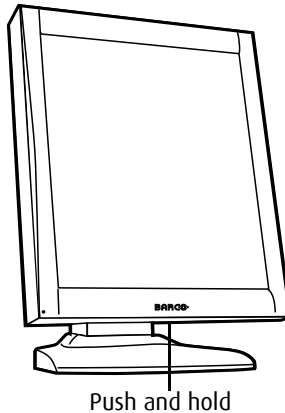



Figure 16: Location of the control wheel

When the display is in stand-by, press the control wheel to switch it back on.

## ON MDNC-2121, MDNC-3121, MDNG-2121 AND MDNG-6121 DISPLAYS:

- 1 When the display is on and no on-screen display is visible, push the control wheel at the front shortly.
- 2 An on-screen message appears, asking to press  to confirm.
- 3 Push the control wheel once more to confirm and switch to stand-by.

- 4 The LED turns orange.

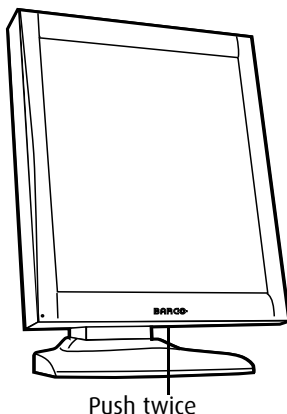


Figure 17: Location of the control wheel

- 5 When the display is in stand-by, press the control wheel shortly to switch it back on.

## Locking and unlocking user controls

The User Controls function allows to disable or enable the control wheel functions. The procedure is slightly different depending on the type of display.

When user controls are disabled, you cannot display and use the on-screen display.

### TO DISABLE USER CONTROLS ON E-3620 AND MDNG-5121:

- 1 Rotate the control wheel to display the on-screen display. The Main Menu appears.
- 2 Rotate the control wheel to select **Settings**.

- 3 Press the control wheel to enter the Settings menu.

SETTINGS	
DPMS	On
Power LED	On
User Controls	On
EXIT	

Figure 18: Settings menu (example)

- 4 Rotate the control wheel to select **User Controls**.
- 5 Press the control wheel to switch from “On” to “Off”.
- 6 Exit the menus.

### TO DISABLE USER CONTROLS ON MDNC-2121, MDNC-3121, MDNG-2121 AND MDNG-6121:

- 1 Rotate the control wheel to display the on-screen display. The Main Menu appears.
- 2 Rotate the control wheel to select **Settings**.
- 3 Press the control wheel to enter the Settings menu.

SETTINGS	
DPMS	On
Power LED	On
User Controls	On
EXIT	

Figure 19: Settings menu (example)

- 4 Rotate the control wheel to select **User Controls**.
- 5 Press the control wheel to activate the User Controls function. This is indicated by the scroll bar that becomes longer.
- 6 Rotate the control wheel to switch from “On” to “Off”.
- 7 Click the control wheel to confirm and return to the Settings menu.
- 8 Exit the menus.



**TO ENABLE USER CONTROLS:**

- 1** Make sure the on-screen display is not on the screen.
- 2** Do not use the control wheel for at least 3 seconds.
- 3** *Only for MDNC-2121, MDNC-3121, MDNG-2121 and MDNG-6121:*  
Press and hold the control wheel for at least 5 seconds.
- 4** Rotate the control wheel 1 step clockwise.
- 5** Press the control wheel 2 times.
- 6** Rotate the control wheel 1 step counterclockwise. .



Steps 4 to 6 must be performed in maximum 3 seconds.

In MDNC-2121, MDNC-3121, MDNG-2121 and MDNG-6121 displays, the user controls are automatically enabled after this procedure. In other Nio displays, you must enable the user controls again in the Settings menu.

# Bringing up the OSD menu

When the OSD is not on the screen, rotate the control wheel to display the OSD.

The main menu appears..


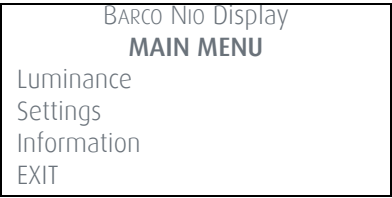


Figure 20: Main menu (example)

## Navigating through the OSD menus

- 1 Rotate the control wheel to move up and down through the menus.
- 2 To go into a submenu, rotate the control wheel to highlight the submenu. Next, click the control wheel to go into the submenu.
- 3 When you click on **Exit**, the OSD exits the submenu and jumps to the previous (higher-level) menu. When you click on **Exit** in the main menu, the OSD is hidden.

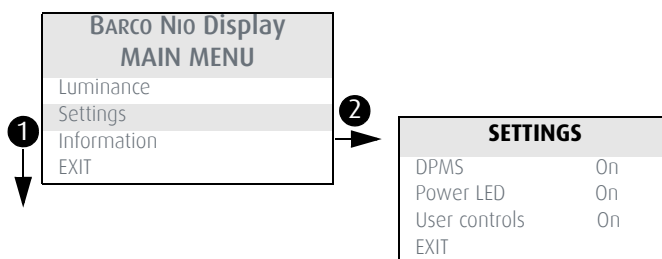

- 1 Rotate the control wheel
- 2 Press the control wheel

Figure 21: OSD menu navigation

## Making selections in the OSD menus

- 1 Navigate through the menus and highlight the parameter you wish to select.
- 2 Click the control wheel to make a selection or to toggle between several settings.

## Changing values in the OSD menus

- 1** Navigate through the menus and highlight the parameter you wish to adjust.
- 2** Click the control wheel to go into edit mode.
- 3** Rotate the wheel to change the value.
- 4** Click the wheel to enter the new value and return to the menu.

# Advanced operation

## On-screen display (OSD)

### Input Selection menu\*

<b>Auto</b>	Detects the incoming video signal automatically and displays it on the screen.
<b>DVI</b>	Always displays the DVI signal on the screen. When no DVI signal available, there will be no image on the screen.
<b>DP</b>	Always displays the DisplayPort signal on the screen. When no DP signal available, there will be no image on the screen.

\* Only in MDNC-2121 displays

# Luminance and color menu



In grayscale Nio displays this menu is called “Luminance”.

<b>Measured luminance**</b>	Indicates the actual luminance measured by the internal sensor. This is a read-only value. It is expressed in a percentage of the calibrated value (100% is calibrated).
<b>Luminance target</b>	Allows to manually adjust the luminance target. See note below.
<b>Color target*</b>	<p>Allows to select from a list of factory-defined and user-defined color targets. You can select from a range of standard Daylight (D) color temperatures, Clearbase, Bluebase and Native.</p> <p>The higher the D-number, the more bluish the color temperature is.</p> <p>‘Native’ is the uncorrected color temperature of the LCD panel.</p>
<b>Viewing mode**</b>	Allows to select a pre-defined viewing mode. <i>Diagnostic</i> is the standard mode for viewing medical images. In <i>Text</i> mode the calibrated luminance is decreased so the viewing conditions are optimal for viewing office documents (text documents, spreadsheets etc.)

\* Only in MDNC-2121 and MDNC-3121 displays.

\*\* Only in MDNC-2121, MDNC-3121, MDNG-2121 and MDNG-6121 displays



When you change the luminance target, the display will adjust its backlight to reach the target. This can be seen in the Measured luminance line.

When the luminance target cannot be reached, e.g., due to aging of the backlight, the Measured luminance line changes to **Minimum value reached** or **Maximum value reached**.

## Display Function menu\*

<b>Display function</b>	<p>Allows to select from a list of pre-defined display functions.</p> <p>If the DICOM DF is selected, an additional function is available.</p>
<b>ALC &amp; DICOM options</b>	<p>Jumps to the ALC &amp; DICOM Options submenu, which allows to edit the settings for the DICOM display function.</p> <p>This function is available only when the DICOM DF is selected.</p>

\*This menu is available only in MDNC-2121, MDNC-3121, MDNG-2121 and MDNG-6121.

## ALC & DICOM OPTIONS SUBMENU\*

<b>Correction value</b>	Shows the ambient light correction value that is taken into account in the calculation of the display function, measured in lux.
<b>DICOM offset</b>	<p>The DICOM DF is recalculated at the moment the DICOM Offset is changed, taking the new DICOM Offset into account.</p> <ul style="list-style-type: none"><li>• When set to “Dark Room”, the ambient light is not taken into account.</li><li>• When set to “Preset”, a preset ambient light value determined by the selected reading room is taken into account.</li></ul>
<b>Reading room</b>	<p>Allows to select from a pre-defined list of reading room types.</p> <p>You must select a reading room that corresponds to the type of room the display is installed in.</p>
<b>Reading room def.</b>	Jumps to the reading room definition submenu, which allows to edit the reading room condition settings.
<b>Calibration info</b>	Jumps to the calculation information submenu, which displays information about the values taken into account to recalculate the DICOM DF.

\*This menu is available only in MDNC-2121, MDNC-3121, MDNG-2121 and MDNG-6121.

**READING ROOM DEFINITION SUBMENU\***

<b>Reading room</b>	Indicates the reading room type you are editing. You can select another room to edit by rotating the control wheel.
<b>Max. ambient light</b>	<p>This indicates the maximum ambient light that corresponds to the selected reading room.</p> <p>If the measured ambient light is higher than the value entered here, you should take measures to darken the room or select another reading room type.</p>
<b>Preset corr. value</b>	<p>This value is taken into account in the calculation of the DICOM DF if DICOM Offset is set to "Preset".</p> <p>To each reading room type corresponds another preset value.</p>

\*This menu is available only in MDNC-2121, MDNC-3121, MDNG-2121 and MDNG-6121.

**CALIBRATION INFORMATION SUBMENU\***

<b>(Preset) Ambient value</b>	Shows the ambient light correction value taken into account to calculate the DICOM DF, expressed in lux.
<b>Bright luminance</b>	Shows the bright luminance value taken into account to calculate the DICOM DF.
<b>Dark luminance</b>	Shows the dark luminance value taken into account to calculate the DICOM DF.
<b>Ambient correction</b>	Shows the ambient light correction value taken into account to calculate the DICOM DF, expressed in $\text{cd/m}^2$ .

\*This menu is available only in MDNC-2121, MDNC-3121, MDNG-2121 and MDNG-6121.



## Settings menu

<b>DPMS</b>	Allows to switch the display power management system on/off. See note below.
<b>Power LED</b>	Allows to switch the power LED's on state on/off. The LED's orange DPMS state is not influenced by this setting. So, when the display goes into power-saving mode, the LED will turn orange, even if it was switched off by this setting
<b>User controls</b>	Allows to disable the control wheel on the front. When switched off, the user cannot display the OSD unit the user controls are enabled again. For more information, <b>Locking and unlocking user controls</b> , on page 23.
<b>Language*</b>	Allows to select another language for the OSD menus.
<b>Auto menu exit**</b>	Allows to switch the automatic menu exit on/off. When switched on, the OSD is closed automatically when left idle for a certain time.
<b>Preferred Input</b>	Select the video resolution you wish the connected graphic board(s) to reproduce. After changing this setting you must reboot the PC. This function takes effect with plug-and-play graphic boards only. In some displays, there is only 1 setting for this function.

\*Only in MDNC-2121, MDNC-3121, MDNG-2121 and MDNG-6121.



The DPMS system will power down the display when the connected computer is left idle for a certain time.



Barco recommends to switch DPMS **on** and to use a good Windows screen saver to prevent image burn-in (image retention) on the LCD panel.

## Information menu

<b>Product</b>	The display type
<b>Serial No</b>	Indicates the display serial number
<b>SW Version</b>	Displays the current internal software version
<b>Display Runtime*</b>	Indicates the total time the display has been operating, including the time in stand-by
<b>Backlight Runtime*</b>	Indicates the total time the display has been operating, excluding the time in stand-by

In some NIO displays these functions are called “Display Lifetime” and “Backlight Lifetime”. Their function, however, is the same as described in the table above.

# Changing Display Functions

## Concepts



Display function selection is possible only in MDNC-3121, MDNG-2121 and MDNG-6121 displays.

- The Display Function menu lets you select a display function (DF).
- When you select the DICOM DF you can define the DICOM DF to be recalculated taking the ambient light into account.
- In this case, the recalculation can be based on an illuminance of 0 lux ("Dark Room") or a preset value ("Preset") determined by the selected reading room.
- In the latter case you must select a reading room that corresponds to the type of room the display is installed in. E.g. if the display is installed in a room where CT scans are observed, the "CT/MR/NM" reading room is the preferred selection.
- The reading room is defined by 2 parameters:
  - **Maximum Ambient Light** defines the maximum light allowed in this type of room. This value can be adjusted within certain limits determined by the selected reading room.
  - **Preset Correction Value** is the predefined correction value for this reading room. This value can be adjusted within certain limits determined by the selected reading room.
- The reading room parameters are pre-defined in the display according to the AAPM (American Association of Physicists in Medicine). However they can be changed within certain limits.

# Display function selection

## TO SELECT A DISPLAY FUNCTION (DF):

- 1 When there is no OSD on the screen, rotate the control wheel to display the OSD. The main menu appears.
- 2 Rotate the control wheel to select the **Display Function** menu. Press the control wheel shortly to enter the menu.
- 3 Rotate the control wheel to select **Display Function**. Press the control wheel shortly to go into edit mode.
- 4 Rotate the control wheel to select the desired DF.
- 5 Press the control wheel shortly to activate the new selection.

You can select the following DFs:

<b>DICOM</b>	Select a DICOM display function for most medical viewing applications.  The DICOM function results in more visible greyscales in the images.
<b>Dynamic Gamma 2.2</b> <b>Dynamic Gamma 1.8</b>	These are gamma functions that are shifted to take into account the non-zero luminance of an LCD panel when driven with a "black" signal.  They are especially useful in CT applications to improve the perception of low Haunsfield values.
<b>Native</b>	If you select Native, the native panel behavior will not be corrected.
<b>Test</b>	This is identical to the Native DF.
<b>Gamma 2.2</b> <b>Gamma 1.8</b>	Select one of these display functions in case the display is to replace a CRT display with a gamma of 1.8 or 2.2 respectively.

User	This is a user-definable DF. Use MediCal QAWeb to define the DF.
------	--

## ALC & DICOM Options

### DICOM offset

#### INTRODUCTION

The DICOM DF will be recalculated *only* when another DICOM Offset is selected. The DICOM DF will be recalculated based on the new offset setting.

#### TO SELECT A DICOM OFFSET:

- 1 Select the DICOM display function. See **To select a display function (DF):**, on page 36.
- 2 Enter the ALC & DICOM submenu.
- 3 In the ALC & DICOM Options submenu, rotate the control wheel to select **DICOM Offset**. Press the control wheel shortly to go into edit mode.
- 4 Rotate the control wheel to select an offset.
- 5 Press the control wheel shortly to activate the selection.

**YOU CAN SELECT THE FOLLOWING OFFSETS:**

<b>Dark room</b>	<p>The DICOM DF will be recalculated taking an ambient light value of 0 lux into account.</p> <p>This means that the measured ambient light does <b>not</b> influence the DICOM DF.</p>
<b>Preset</b>	<p>The DICOM DF will be recalculated taking a preset ambient light value into account. This value is determined by the selected Reading Room and can be adjusted manually within certain limits.</p>

**Reading room selection**

This function is available only when the DICOM display function is selected.

**TO SELECT A READING ROOM SETTING:**

- 1** Select the DICOM display function. See **To select a display function (DF):**, on page 36.
- 2** Enter the ALC & DICOM Options submenu.
- 3** In the ALC & DICOM Options submenu, rotate the control wheel to select **Reading Room**. Press the control wheel shortly to go into edit mode.
- 4** Rotate the control wheel to select the reading room that corresponds to the room type the display is installed in. Return to the menu.
- 5** You can select the following reading room types:

<b>CR/DR/MAMMO</b>	<p>Corresponds to light conditions in diagnostic reading rooms for computed radiology, digital radiology or mammography. This setting has the lowest maximum ambient light.</p>
--------------------	---

<b>CT/MR/NM</b>	Corresponds to light conditions in diagnostic reading rooms for computed tomography, magnetic resonance or nuclear medicine scans.
<b>Office</b>	Corresponds to light conditions in office rooms.
<b>Clinical</b>	Corresponds to light conditions in diagnostic reading rooms for clinical viewing.
<b>Emergency</b>	Corresponds to light conditions in emergency rooms.
<b>Operating</b>	Corresponds to light conditions in operating rooms. This setting has the highest maximum ambient light.

## Reading room definition

This function is available only when the DICOM display function is selected.

### TO MODIFY A READING ROOM DEFINITION:

- 1 Select the DICOM display function. See **To select a display function (DF):**, on page 36.
- 2 Enter the ALC & DICOM Options submenu.
- 3 In the ALC & DICOM Options submenu, rotate the control wheel to select **Reading Room Def**. Press the control wheel shortly to open the submenu.
- 4 Select the reading room definition you wish to modify. Thereto, rotate the control wheel to select **Reading Room**. Press the control wheel shortly to go into edit mode. Rotate the control wheel to select the reading room. Press the control wheel shortly to confirm the selection.

- 5 To modify the *maximum ambient light* for this reading room, rotate the control wheel to select **Max. Ambient Light**. Press the control wheel shortly to go into edit mode. Rotate the control wheel to change the value. Press the control wheel shortly to save the changes.

You can change this value within certain limits determined by the selected reading room.

- 6 To modify the *Preset correction value* for this reading room, rotate the control wheel to select **Preset Corr. Value**. Press the control wheel shortly to go into edit mode. Rotate the control wheel to change the value. Press the control wheel shortly to save the changes.

You can change this value within certain limits determined by the selected reading room.

## Switching DPMS on/off

### TO SWITCH DPMS ON/OFF:

- 1 When there is no OSD on the screen, rotate the control wheel. The main OSD menu appears.
- 2 Rotate the control wheel to select the **Settings** menu. Press the control wheel shortly to open the menu.
- 3 Rotate the control wheel to select **DPMS**.
- 4 For E-3620 and MDNG-5121: Press the control wheel to change the DPMS setting.  
For MDNC-3121, MDNG-2121 and MDNG-6121:  
Press the control wheel shortly to go into edit mode. Next, rotate the control wheel to change the setting. At last press the control wheel again to save the changes.



# Cleaning your display

## Front glass

Clean the front glass using a sponge, cleaning cloth or soft tissue, **lightly** moistened with one of the following tested products:

- Flux
- Windex Glass Plus
- Bohle glass cleaner
- Mr. Proper
- Pril
- Ajax glass cleaner
- Sidolin glass cleaner
- 25% ethanol solution (EtOH)

In case none of the above cleaning products is available, use plain water.



Take care not to damage or scratch the front glass. Be careful with rings or other jewelry and do not apply excessive pressure on the front glass.



Do not apply or spray liquid directly to the front glass as excess liquid may cause damage to internal electronics. Instead, apply the liquid to the cleaning cloth.



Do not use following products:

- Alcohol/solvents at higher concentration > 5%
- Strong alkalis lye, strong solvents
- Acid
- Detergents with fluoride
- Detergents with ammonia
- Detergents with abrasives
- Steel wool
- Sponge with abrasives
- Steel blades
- Cloth with steel thread

## Cabinet

- 1 Clean the cabinet using a soft cotton cloth, **lightly** moistened with a recognized cleaning product for medical equipment.

The cabinet has been tested for resistance to the following products:

- Cidex, Betadine
- Alcohol (Isopropyl and Ethyl)
- Ammonia-based cleaners (Windex)
- Aquasonic Gel

In case none of the above cleaning products is available, use plain water.



Do not apply or spray liquid directly to the cabinet as excess liquid may cause damage to internal electronics. Instead, apply the liquid to the cleaning cloth.

- 2 Repeat with water only.
- 3 Wipe dry with a dry cloth.

# Important information

## Safety information

### GENERAL RECOMMENDATIONS


Read the safety and operating instructions before operating the equipment.

Retain safety and operating instructions for future reference.

Adhere to all warnings on the equipment and in the operating instructions manual.

Follow all instructions for operation and use.

### ELECTRICAL SHOCK

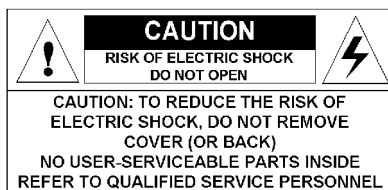

### TYPE OF PROTECTION (ELECTRICAL):

Equipment with external power supply: Class I equipment

### DEGREE OF SAFETY (FLAMMABLE ANESTHETIC MIXTURE):

Equipment not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide.

### NON-PATIENT CARE EQUIPMENT

Equipment primarily for use in a health care facility that is intended for use where contact with a patient is unlikely.

## **POWER CONNECTION - EQUIPMENT WITH EXTERNAL 12 VDC POWER SUPPLY**

- Power requirements: The equipment must be powered using the delivered medical approved 12 VDC power supply.
- The medical approved DC power supply must be powered by the AC mains voltage (protective earth terminal).

## **POWER CONNECTION - EQUIPMENT WITH EXTERNAL 24 VDC POWER SUPPLY**

- Power requirements: The equipment must be powered using the delivered medical approved 24 VDC SELV power supply.
- The medical approved DC power supply must be powered by the AC mains voltage (protective earth terminal).

## **TRANSIENT OVER-VOLTAGE**

If the equipment is not used for a long time, disconnect it from the AC inlet to avoid damage by transient over-voltage.

To fully disengage the power to the unit, please disconnect the power cord from the AC inlet.

## **POWER CORDS:**

- Utilize a UL-listed detachable power cord, 3-wire, type SJ or equivalent, 18 AWG min., rated 300 V min., provided with a hospital-grade type plug 5-15P configuration for 120V application, or 6-15P for 240V application.
- Do not overload wall outlets and extension cords as this may result in fire or electric shock.
- Mains lead protection (U.S.: Power cord): Power cords should be routed so that they are not likely to be walked upon or pinched by items placed upon or against them, paying particular attention to cords at plugs and receptacles.

## **WATER AND MOISTURE**

Never expose the equipment to rain or moisture.

Never use the equipment near water - e.g. near a bathtub, washbasin, swimming pool, kitchen sink, laundry tub or in a wet basement.

## **VENTILATION**

Do not cover or block the ventilation openings in the cover of the set. When installing the equipment in a cupboard or another closed location, heed the necessary space between the set and the sides of the cupboard.

## **INSTALLATION**

Place the equipment on a flat, solid and stable surface that can support the weight of at least 3 equipments. If you use an unstable cart or stand, the equipment may fall, causing serious injury to a child or adult, and serious damage to the equipment.

### **THIS APPARATUS CONFORMS TO:**

CE0120 (MDD 93/42/EEC class IIb product), IEC 60601-1, UL 60601-1,

CAN/CSA C22.2 No. 601.01-M90 (c-UL),

CCC GB4943-1995 (IEC 60950-1).

**MDNC-2121:** Medical equipment with respect to electric shock, fire and mechanical hazards only in accordance with UL 60601-1. and CAN/CSA C22.2 NO. 601.1.

### **NATIONAL SCANDINAVIAN DEVIATIONS FOR CL. 1.7.2:**

Finland: "Laite on liitettävä suojamaadoituskoskettimilla varustettuun pistorasiaan"

Norway: "Apparatet må tilkoples jordet stikkontakt"

Sweden: "Apparaten skall anslutas till jordat uttag"

# Environmental information

## DISPOSAL INFORMATION (WASTE ELECTRICAL AND ELECTRONIC EQUIPMENT)



This symbol on the product indicates that, under European Directive 2002/96/EC governing waste from electrical and electronic equipment, this product must not be disposed of with other municipal waste. Please dispose of your waste equipment by handing it over to a designated collection point for the recycling of waste electrical and electronic equipment. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate these items from other types of waste and recycle them responsibly to promote the sustainable reuse of material resources.

For more information about recycling of this product, please contact your local city office, your municipal waste disposal service or the shop where you purchased the product.

## MERCURY NOTICE

This Barco product consists of materials that may contain mercury, which must be recycled or disposed of in accordance with local, state, or country laws:

- Within this system, the backlight lamps in the monitor display contain mercury



## TURKEY RoHS COMPLIANCE



Republic of Turkey: In conformity with the EEE Regulation  
Türkiye Cumhuriyeti: EEE Yönetmeliğine Uygundur

## 中国大陆 RoHS (CHINESE MAINLAND RoHS)

根据中国大陆《电子信息产品污染控制管理办法》（也称为中国大陆 RoHS），以下部分列出了 Barco 产品中可能包含的有毒和 / 或有害物质的名称和含量。中国大陆 RoHS 指令包含在中国信息产业部 MCV 标准：“电子信息产品中有毒物质的限量要求”中。

According to the "China Administration on Control of Pollution Caused by Electronic Information Products" (Also called RoHS of Chinese Mainland), the table below lists the names and contents of toxic and/or hazardous substances that Barco's product may contain. The RoHS of Chinese Mainland is included in the MCV standard of the Ministry of Information Industry of China, in the section "Limit Requirements of toxic substances in Electronic Information Products".

零件项目 (名称) Component Name	有毒有害物质或元素 Hazardous Substances or Elements					
	铅 (Pb)	汞 (Hg)	镉 (Cd)	六价铬 (Cr6+)	多溴联苯 (PBB)	多溴二苯醚 (PBDE)
印制电路配件 Printed Circuit Assemblies	X	0	0	0	0	0
液晶面板 LCD panel	X	X	0	0	0	0
外接电 (线) 缆 External cables	X	0	0	0	0	0
内部线路 Internal wiring	0	0	0	0	0	0
金属外壳 Metal enclosure	0	0	0	0	0	0
塑胶外壳 Plastic enclosure	0	0	0	0	0	0
散热片 (器) Heatsinks	0	0	0	0	0	0
风扇 Fan	0	0	0	0	0	0

电源供应器 Power supply unit	X	0	0	0	0	0
文件说明书 Paper manuals	0	0	0	0	0	0
光盘说明书 CD manual	0	0	0	0	0	0
<p>O: 表示该有毒有害物质在该部件所有均质材料中的含量均在 SJ/T 11363-2006 标准规定的限量要求以下。 O: Indicates that this toxic or hazardous substance contained in all of the homogeneous materials for this part is below the limit requirement in SJ/T11363-2006.</p> <p>X: 表示该有毒有害物质至少在该部件的某一均质材料中的含量超出 SJ/T 11363-2006 标准规定的限量要求。 X: Indicates that this toxic or hazardous substance contained in at least one of the homogeneous materials used for this part is above the limit requirement in SJ/T11363 2006.</p>						

在中国大陆销售的相应电子信息产品（EIP）都必须遵照中国大陆《电子信息产品污染控制标识要求》标准贴上环保使用期限（EFUP）标签。Barco 产品所采用的 EFUP 标签（请参阅实例，徽标内部的编号使用于制定产品）基于中国大陆的《电子信息产品环保使用期限通则》标准。

All Electronic Information Products (EIP) that are sold within Chinese Mainland must comply with the "Electronic Information Products Pollution Control Labeling Standard" of Chinese Mainland, marked with the Environmental Friendly Use Period (EFUP) logo. The number inside the EFUP logo that Barco uses (please refer to the photo) is based on the "Standard of Electronic Information Products Environmental Friendly Use Period" of Chinese Mainland.





# Regulatory compliance information

## **INDICATIONS FOR USE: E-3620, MDNG-2121, MDNC-3121, MDNC-2121**

The products are intended to be used in displaying and viewing digital images, for review and analysis by trained medical practitioners. These devices must not be used in primary image diagnosis in mammography.

Caution (USA): Federal law restricts this device to sale by or no the order of a physician or a practitioner trained on its use.

## **INDICATIONS FOR USE: MDNG-5121, MDNG-6121**

The products are intended to be used in displaying and viewing digital images, including digital mammography, for review and analysis by trained medical practitioners.

Caution (USA): Federal law restricts this device to sale by or no the order of a physician or a practitioner trained on its use.

## **FCC CLASS B**

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:










- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.







### CANADIAN NOTICE

This ISM device complies with Canadian ICES-001.  
 Cet appareil ISM est conforme à la norme NMB-001 du Canada.

## Explanation of symbols

On the display or power supply, you may find the following symbols (nonrestrictive list):

	Indicates compliance to the essential requirements of the Directive 93/42/EEC
	Indicates compliance with Part 15 of the FCC rules (Class A or Class B)
 or 	Indicates the display is approved according to the UL regulations
 or 	Indicates the display is approved according to the c-UL regulations
	Indicates the display is approved according to the DEMKO regulations
	Indicates the display is approved according to the CCC regulations
	Indicates the display is approved according to the VCCI regulations

	Indicates the USB connectors on the display
	Indicates the manufacturing date
	Indicates the temperature limitations for the display to operate within specs
	Indicates the display serial no.
	Consult the operating instructions
	Indicates this apparatus must not be thrown in the trash but must be recycled, according to the European WEEE (Waste Electrical and Electronic Equipment) directive

# Legal disclaimer

## DISCLAIMER NOTICE

Although every attempt has been made to achieve technical accuracy in this document, we assume no responsibility for errors that may be found. Our goal is to provide you with the most accurate and usable documentation possible; if you discover errors, please let us know.

Barco software products are the property of Barco. They are distributed under copyright by Barco N.V. or Barco, Inc., for use only under the specific terms of a software license agreement between Barco N.V. or Barco Inc. and the licensee. No other use, duplication, or disclosure of a Barco software product, in any form, is authorized.

The specifications of Barco products are subject to change without notice.

## TRADEMARKS

All trademarks and registered trademarks are property of their respective owners.

## COPYRIGHT NOTICE

This document is copyrighted. All rights are reserved. Neither this document, nor any part of it, may be reproduced or copied in any form or by any means - graphical, electronic, or mechanical including photocopying, taping or information storage and retrieval systems - without written permission of Barco.

© 2012 Barco N.V. All rights reserved.

# Technical specifications

## Nio 2MP

Product acronym	MDNG-2121
Screen technology	TFT AM LCD Dual Domain IPS
Active screen size (diagonal)	540 mm (21.3")
Active screen size (H x V)	432 x 324 mm (17.0 x 12.8")
Aspect ratio (H:V)	4:3
Resolution	2MP (1600 x 1200)
Pixel pitch	0.2700 mm
Color imaging	No
Gray imaging	Yes
Number of grayscales (LUT in/LUT out)	1024 grey levels (10/12)
Viewing angle (H, V)	170°
Uniform Luminance Technology (ULT)	Yes
Per Pixel Uniformity (PPU)	No
Ambient Light Compensation (ALC)	No
Backlight Output Stabilization (BLOS)	Yes
I-Guard	No
Maximum luminance*	1650 cd/m²
DICOM calibrated luminance (ULT off)*	500 cd/m²
Contrast ratio (typical)*	850:1

<b>Response time (Tr + Tf)</b>	23 ms
<b>Scanning frequency (H; V)</b>	15-129 kHz; 24-100 Hz
<b>Housing color</b>	Gray
<b>Video input signals</b>	DVI-D Dual Link
<b>USB ports</b>	1 upstream (endpoint), 2 downstream
<b>USB standard</b>	2.0
<b>Power requirements (nominal)</b>	100-240V
<b>Power consumption (nominal)</b>	65W
<b>Power save mode</b>	Yes
<b>Power management</b>	DVI-DMPM
<b>Dot clock</b>	280 MHz
<b>OSD languages</b>	English, French, German, Spanish, Italian
<b>Dimensions with stand (W x H x D)</b>	Portrait: 385 x 525~585 x 250 mm Landscape: 485 x 475~535 x 250 mm
<b>Dimensions w/o stand (W x H x D)</b>	Portrait: 385 x 485 x 115 mm Landscape: 485 x 385 x 115 mm
<b>Dimensions packaged (W x H x D)</b>	500 x 675 x 320 mm
<b>Net weight with stand</b>	12.2 kg
<b>Net weight w/o stand</b>	8.14 kg
<b>Net weight packaged with stand</b>	17.19 kg
<b>Net weight packaged w/o stand</b>	13.15 kg
<b>Height adjustment range</b>	60 mm
<b>Tilt</b>	-5° / +30°

<b>Swivel</b>	-45° / +45°
<b>Pivot</b>	90°
<b>Mounting standard</b>	VESA (100 mm)
<b>Screen protection</b>	Protective, non-reflective PMMA cover
<b>Recommended modalities</b>	CT, MR, US, DR, CR, NM, Film
<b>Certifications</b>	BSMI, CCC - GB4943-2001 - GB9254-1998 - GB9254-2008 - GB17625.1-2003, DEMKO EN60601-1:1990 + A1:1993 + A2:1995, IEC 60601-1:1998 + A1:91 + A2:95, UL 60601-1 1st edition, CAN/CSA-CS22.2 No.601.1-M90, FDA510K, RoHs, MDD, WEEE,CE 0120, VCCI, FCC Class B, KETI
<b>Supplied accessories</b>	User Guide Quick-Installation Sheet Video cable (DVI Dual Link) Mains cables (UK, European (CEBEC/KEMA), USA (UL/CSA; adaptor plug NEMA 5-15P), Chinese (CCC)) USB 2.0 cable External power supply
<b>Optional accessories</b>	NA
<b>QA software</b>	MediCal QAWeb
<b>Unit per pallet</b>	NA
<b>Pallet dimensions</b>	NA
<b>Warranty</b>	5 years
<b>Operating temperature</b>	0°C to 35°C (15°C to 35°C within specs)
<b>Storage temperature</b>	-20°C to 60°C
<b>Operating humidity</b>	8% - 80% (non-condensing)
<b>Storage humidity</b>	5% - 95% (non-condensing)
<b>Operating altitude</b>	3000 m
<b>Storage altitude</b>	7500 m

\* In dark reading room conditions (0 lux)

# Nio 3MP

Product acronym	E-3620 (MA)
Screen technology	TFT AM LCD Dual Domain IPS
Active screen size (diagonal)	528 mm (20.8")
Active screen size (H x V)	424 x 318 mm (16.7 x 12.5")
Aspect ratio (H:V)	4:3
Resolution	3MP (2048 x 1536)
Pixel pitch	0.207 mm
Color imaging	No
Gray imaging	Yes
Number of grayscales (LUT in/LUT out)	1024 grey levels (10/10)
Viewing angle (H, V)	170°
Uniform Luminance Technology (ULT)	Yes
Per Pixel Uniformity (PPU)	No
Ambient Light Compensation (ALC)	No
Backlight Output Stabilization (BLOS)	Yes
I-Guard	No
Maximum luminance*	1000 cd/m <sup>2</sup>
DICOM calibrated luminance (ULT off)*	500 cd/m <sup>2</sup>
Contrast ratio (typical)*	900:1
Response time (Tr + Tf)	50 ms
Scanning frequency (H; V)	30-124 kHz; 50-85 Hz



<b>Housing color</b>	Gray
<b>Video input signals</b>	DVI-D Dual Link
<b>USB ports</b>	1 upstream, 2 downstream
<b>USB standard</b>	1.1
<b>Power requirements (nominal)</b>	100-240V
<b>Power consumption (nominal)</b>	54W
<b>Power save mode</b>	Yes
<b>Power management</b>	DVI-DMPM
<b>Dot clock</b>	280 MHz
<b>OSD languages</b>	English
<b>Dimensions with stand (W x H x D)</b>	Portrait: 385 x 525~585 x 250 mm Landscape: 485 x 475~535 x 250 mm
<b>Dimensions w/o stand (W x H x D)</b>	Portrait: 385 x 485 x 115 mm Landscape: 485 x 385 x 115 mm
<b>Dimensions packaged (W x H x D)</b>	500 x 675 x 320 mm
<b>Net weight with stand</b>	13 kg
<b>Net weight w/o stand</b>	9 kg
<b>Net weight packaged with stand</b>	15 kg
<b>Net weight packaged w/o stand</b>	11 kg
<b>Height adjustment range</b>	60 mm
<b>Tilt</b>	-5° / +30°
<b>Swivel</b>	-45° / +45°
<b>Pivot</b>	90°

<b>Mounting standard</b>	VESA (100 mm)
<b>Screen protection</b>	Protective, non-reflective PMMA cover
<b>Recommended modalities</b>	CT, MR, US, DR, CR, NM, Film
<b>Certifications</b>	CCC - GB4943-2001 - GB9254-1998 (Class A) - GB9254-2008 (Class A) - GB17625.1-2003, DEMKO EN60601-1:1990 + A1:1993 + A2:1995 + A11: 1993 + A12: 1993 + A13:1996, IEC 60601-1:1988 + A1:1991 + A2:1995 2nd edition, UL 60601-1 1st edition 2006-04-26, CAN/CSA-C22.2 No.601.1-M90 2005, FDA510K, RoHS, MDD, WEEE, CE 0120, VCCI, FCC Class B, KETI
<b>Supplied accessories</b>	User Guide Quick-Installation Sheet Video cable (DVI Dual Link) Mains cables (UK, European (CEBEC/KEMA), USA (UL/CSA; adaptor plug NEMA 5-15P), Chinese (CCC)) External power supply
<b>Optional accessories</b>	NA
<b>QA software</b>	MediCal QAWeb
<b>Unit per pallet</b>	NA
<b>Pallet dimensions</b>	NA
<b>Warranty</b>	5 years
<b>Operating temperature</b>	0°C to 40°C (15°C to 35°C within specs)
<b>Storage temperature</b>	-20°C to 60°C
<b>Operating humidity</b>	8% - 80% (non-condensing)
<b>Storage humidity</b>	5% - 95% (non-condensing)
<b>Operating altitude</b>	3000 m
<b>Storage altitude</b>	7500 m

\* In dark reading room conditions (0 lux)

# Nio 5MP

Product acronym	MDNG-5121
Screen technology	TFT AM LCD Dual Domain IPS
Active screen size (diagonal)	541 mm (21.3")
Active screen size (H x V)	422 x 338 mm (16.5 x 13.3")
Aspect ratio (H:V)	5:4
Resolution	5MP (2560 x 2048)
Pixel pitch	0.165 mm
Color imaging	No
Gray imaging	Yes
Number of grayscales (LUT in/LUT out)	10 bit/10 bit
Viewing angle (H, V)	170°
Uniform Luminance Technology (ULT)	No
Per Pixel Uniformity (PPU)	No
Ambient Light Compensation (ALC)	No
Backlight Output Stabilization (BLOS)	Yes
I-Guard	No
Maximum luminance*	700 cd/m²
DICOM calibrated luminance (ULT off)*	500 cd/m²
Contrast ratio (typical)*	800:1
Response time (Tr + Tf)	50 ms
Scanning frequency (H; V)	30-124 kHz; 45-85 Hz

<b>Housing color</b>	Gray
<b>Video input signals</b>	DVI-D Dual Link
<b>USB ports</b>	1 upstream (endpoint), 2 downstream
<b>USB standard</b>	1.1
<b>Power requirements (nominal)</b>	100-250V
<b>Power consumption (nominal)</b>	61W
<b>Power save mode</b>	Yes
<b>Power management</b>	DVI-DMPM
<b>Dot clock</b>	330 MHz
<b>OSD languages</b>	English
<b>Dimensions with stand (W x H x D)</b>	Portrait: 408 x 489~549 x 250 mm Landscape: 492 x 531~591 x 250 mm
<b>Dimensions w/o stand (W x H x D)</b>	Portrait: 408 x 492 x 115 mm Landscape: 492 x 408 x 115 mm
<b>Dimensions packaged (W x H x D)</b>	500 x 320 x 675 mm
<b>Net weight with stand</b>	12.5 kg
<b>Net weight w/o stand</b>	8.5 kg
<b>Net weight packaged with stand</b>	15.2 kg
<b>Net weight packaged w/o stand</b>	NA
<b>Height adjustment range</b>	60 mm
<b>Tilt</b>	-5° / +30°
<b>Swivel</b>	-45° / +45°
<b>Pivot</b>	90°

<b>Mounting standard</b>	VESA (100 mm)
<b>Screen protection</b>	Protective, non-reflective PMMA cover
<b>Recommended modalities</b>	CT, MR, US, DR, CR, NM, Film
<b>Certifications</b>	BSMI, CCC - GB4943-2001 - GB9254-1998 (Class A) - GB9254-2008 (Class A) - GB17625.1-2003, DEMKO EN60601-1:1990 + A1:1993 + A2:1995, IEC 60601-1:1998 + A1:1991 + A2:1995 2nd edition, IEC 60601-1 (ed.2) am1 & am2, UL 60601-1 1st edition 2006-04-26, CAN/CSA-C22.2 No.601.1-M90 2005, FDA510K, RoHS, MDD, WEEE, CE 0120, VCCI, FCC Class B, KETI
<b>Supplied accessories</b>	User Guide Quick-Installation Sheet Video cable (DVI Dual Link) Mains cables (UK, European (CEBEC/KEMA), USA (UL/CSA; adaptor plug NEMA 5-15P), Chinese (CCC)) USB 2.0 cable External power supply
<b>Optional accessories</b>	NA
<b>QA software</b>	MediCal QAWeb
<b>Unit per pallet</b>	NA
<b>Pallet dimensions</b>	NA
<b>Warranty</b>	5 years
<b>Operating temperature</b>	0°C to 40°C (15°C to 35°C within specs)
<b>Storage temperature</b>	-20°C to 60°C
<b>Operating humidity</b>	8% - 80% (non-condensing)
<b>Storage humidity</b>	5% - 95% (non-condensing)
<b>Operating altitude</b>	3000 m
<b>Storage altitude</b>	7500 m

\* In dark reading room conditions (0 lux)

# Nio COLOR 3MP HB

Product acronym	MDNC-3121
Screen technology	TFT AM Color LCD IPS
Active screen size (diagonal)	540 mm (21.3")
Active screen size (H x V)	433 x 325 mm (17.0 x 12.8")
Aspect ratio (H:V)	4:3
Resolution	3MP (2048 x 1536)
Pixel pitch	0.2115 mm
Color imaging	Yes
Gray imaging	No
Color support	24 bit
Viewing angle (H, V)	176°
Uniform Luminance Technology (ULT)	Yes
Per Pixel Uniformity (PPU)	No
Ambient Light Compensation (ALC)	No
Backlight Output Stabilization (BLOS)	Yes
I-Guard	No
Maximum luminance*	800 cd/m²
DICOM calibrated luminance (ULT off)*	400 cd/m²
Contrast ratio (typical)*	750:1
Response time (Tr + Tf)	24 ms
Scanning frequency (H; V)	15-128 kHz; 25-98 Hz
Housing color	Gray

<b>Video input signals</b>	DVI-D Dual Link
<b>USB ports</b>	1 upstream (endpoint), 2 downstream
<b>USB standard</b>	2.0
<b>Power requirements (nominal)</b>	100-250V
<b>Power consumption (nominal)</b>	72W
<b>Power save mode</b>	Yes
<b>Power management</b>	DVI-DMPM
<b>Dot clock</b>	165 MHz
<b>OSD languages</b>	English, French, German, Spanish, Italian
<b>Dimensions with stand (W x H x D)</b>	Portrait: 385 x 525~585 x 250 mm Landscape: 485 x 475~535 x 250 mm
<b>Dimensions w/o stand (W x H x D)</b>	Portrait: 385 x 485 x 115 mm Landscape: 485 x 385 x 115 mm
<b>Dimensions packaged (W x H x D)</b>	500 x 320 x 675 mm
<b>Net weight with stand</b>	13 kg
<b>Net weight w/o stand</b>	9 kg
<b>Net weight packaged with stand</b>	15 kg
<b>Net weight packaged w/o stand</b>	NA
<b>Height adjustment range</b>	60 mm
<b>Tilt</b>	-5° / +30°
<b>Swivel</b>	-45° / +45°
<b>Pivot</b>	90°
<b>Mounting standard</b>	VESA (100 mm)

<b>Screen protection</b>	Protective, non-reflective PMMA cover
<b>Recommended modalities</b>	CT, MR, US, DR, CR, NM, Film
<b>Certifications</b>	B5MI, CCC - GB4943-2001 - GB9254-2008 - GB17625.1-2003, DEMKO EN60601-1:1990 + A1:1993 + A2:1995 + A13:1996, IEC 60601-1:1988 + A1:1991 + A2:1995 2nd edition, UL 60601-1 1st edition 2006-04-26, CAN/CSA-C22.2 No.601.1-M90 2005, FDA510K, RoHs, MDD, WEEE, CE 0120, VCCI, FCC Class B, KETI
<b>Supplied accessories</b>	User Guide Quick-Installation Sheet Video cable (DVI Dual Link) Mains cables (UK, European (CEBEC/KEMA), USA (UL/CSA; adaptor plug NEMA 5-15P), Chinese (CCC)) USB 2.0 cable External power supply
<b>Optional accessories</b>	NA
<b>QA software</b>	MediCal QAWeb
<b>Unit per pallet</b>	NA
<b>Pallet dimensions</b>	NA
<b>Warranty</b>	5 years
<b>Operating temperature</b>	0°C to 40°C (15°C to 35°C within specs)
<b>Storage temperature</b>	-20°C to 60°C
<b>Operating humidity</b>	8% - 70% (non-condensing)
<b>Storage humidity</b>	5% - 95% (non-condensing)
<b>Operating altitude</b>	3600 m
<b>Storage altitude</b>	7500 m

\* In dark reading room conditions (0 lux)



# Nio 5MP

Product acronym	MDNG-6121
Screen technology	TFT AM LCD Dual Domain IPS
Active screen size (diagonal)	541 mm (21.3")
Active screen size (H x V)	432,6 x 324,5 mm (17.0 x 12.8")
Aspect ratio (H:V)	4:3
Resolution	5MP (2800 x 2100)
Pixel pitch	0.1540 mm
Color imaging	No
Gray imaging	Yes
Number of grayscales (LUT in/LUT out)	1024 grey levels (10/12)
Viewing angle (H, V)	170°
Uniform Luminance Technology (ULT)	Yes
Per Pixel Uniformity (PPU)	No
Ambient Light Compensation (ALC)	No
Backlight Output Stabilization (BLOS)	Yes
I-Guard	No
Maximum luminance*	1100 cd/m <sup>2</sup>
DICOM calibrated luminance (ULT off)*	500 cd/m <sup>2</sup>
Contrast ratio (typical)*	1000:1
Response time (Tr + Tf)	32 ms
Scanning frequency (H; V)	15-129 kHz; 24-100 Hz

<b>Housing color</b>	Gray
<b>Video input signals</b>	DVI-D Dual Link
<b>USB ports</b>	1 upstream (endpoint), 2 downstream
<b>USB standard</b>	2.0
<b>Power requirements (nominal)</b>	100-240V
<b>Power consumption (nominal)</b>	78W
<b>Power save mode</b>	Yes
<b>Power management</b>	DVI-DMPM
<b>Dot clock</b>	280 MHz
<b>OSD languages</b>	English, French, German, Spanish, Italian
<b>Dimensions with stand (W x H x D)</b>	Portrait: 385 x 525~585 x 250 mm Landscape: 485 x 475~535 x 250 mm
<b>Dimensions w/o stand (W x H x D)</b>	Portrait: 385 x 485 x 115 mm Landscape: 485 x 385 x 115 mm
<b>Dimensions packaged (W x H x D)</b>	500 x 320 x 675 mm
<b>Net weight with stand</b>	12.3 kg
<b>Net weight w/o stand</b>	8.3 kg
<b>Net weight packaged with stand</b>	17.3 kg
<b>Net weight packaged w/o stand</b>	13.3 kg
<b>Height adjustment range</b>	60 mm
<b>Tilt</b>	-5° / +30°
<b>Swivel</b>	-45° / +45°
<b>Pivot</b>	90°


<b>Mounting standard</b>	VESA (100 mm)
<b>Screen protection</b>	Protective, non-reflective PMMA cover
<b>Recommended modalities</b>	CT, MR, US, DR, CR, NM, Film
<b>Certifications</b>	BSMI, CCC - GB4943-2001 - GB9254-1998 - GB9254-2008 - GB17625.1-2003, DEMKO EN60601-1:1990 + A1:1993 + A2:1995 + A13:1996, IEC 60601-1:1988 + A1:1991 + A2:1995 2nd edition, UL 60601-1 1st edition 2006-04-26, CAN/CSA-C22.2 No.601.1-M90 2005, FDA510K, RoHs, MDD, WEEE, CE 0120, VCCI, FCC Class B, KETI
<b>Supplied accessories</b>	User Guide Quick-Installation Sheet Video cable (DVI Dual Link) Mains cables (UK, European (CEBEC/KEMA), USA (UL/CSA; adaptor plug NEMA 5-15P), Chinese (CCC)) USB 2.0 cable External power supply
<b>Optional accessories</b>	NA
<b>QA software</b>	MediCal QAWeb
<b>Unit per pallet</b>	NA
<b>Pallet dimensions</b>	NA
<b>Warranty</b>	5 years
<b>Operating temperature</b>	0°C to 35°C (15°C to 30°C within specs)
<b>Storage temperature</b>	-20°C to 60°C
<b>Operating humidity</b>	8% - 70% (non-condensing)
<b>Storage humidity</b>	5% - 95% (non-condensing)
<b>Operating altitude</b>	3000 m
<b>Storage altitude</b>	7500 m

\* In dark reading room conditions (0 lux)

# NIO COLOR 2MP HB

Product acronym	MDNC-2121
Screen technology	TFT AM Color LCD IPS
Active screen size (diagonal)	541 mm (21.3")
Active screen size (H x V)	432 x 324 mm (17.0 x 12.8")
Aspect ratio (H:V)	4:3
Resolution	2MP (1600 x 1200)
Pixel pitch	0.2700 mm (0.01063")
Color imaging	Yes
Gray imaging	No
Color support	24 bit
Viewing angle (H, V)	176°
Uniform Luminance Technology (ULT)	Yes
Per Pixel Uniformity (PPU)	No
Ambient Light Compensation (ALC)	No
Backlight Output Stabilization (BLOS)	Yes
I-Guard	No
Maximum luminance*	750 cd/m²
DICOM calibrated luminance (ULT off)*	400 cd/m²
Contrast ratio (typical)*	1050:1
Response time (Tr + Tf)	35 ms
Scanning frequency (H; V)	15-129 kHz; 24-100 Hz
Housing color	Gray

<b>Video input signals</b>	DVI-D Single Link (1x), DisplayPort (1x)
<b>USB ports</b>	1 upstream (endpoint), 2 downstream
<b>USB standard</b>	2.0
<b>Power source</b>	Input for 24 VDC power supply unit: 90-264 VAC Input for display: 24 VDC For use only power: Sinpro HPU100-108
<b>Power consumption (nominal)</b>	66W
<b>Power save mode</b>	Yes
<b>Power management</b>	DVI-DMPM
<b>Dot clock</b>	165 MHz
<b>OSD languages</b>	English, French, German, Spanish, Italian
<b>Dimensions with stand (W x H x D)</b>	Portrait: 385 x 525~585 x 250 mm Landscape: 485 x 475~535 x 250 mm
<b>Dimensions w/o stand (W x H x D)</b>	Portrait: 385 x 485 x 115 mm Landscape: 485 x 385 x 115 mm
<b>Dimensions packaged (W x H x D)</b>	500 x 320 x 675 mm
<b>Net weight with stand</b>	13 kg
<b>Net weight w/o stand</b>	9 kg
<b>Net weight packaged with stand</b>	15 kg
<b>Net weight packaged w/o stand</b>	NA
<b>Height adjustment range</b>	60 mm
<b>Tilt</b>	-5° / +30°
<b>Swivel</b>	-45° / +45°
<b>Pivot</b>	90°
<b>Mounting standard</b>	VESA (100 mm)

<b>Screen protection</b>	Protective, non-reflective PMMA cover
<b>Recommended modalities</b>	CT, MR, US, DR, CR, NM, Film
<b>Certifications</b>	UL60601-1, CB60601-1, CSA C22.2 No 601.1, EN60601-1-2, EN55011 Limit B, ICES-001 B, EN61000-3-2/3, EN61000-4-2/3/4/5/6/8/11, VCCI, CE0120, KETI, CCC, BSMI, Demko, FCCB, FDA 510K, CMDCAS, MDD, WEEE, RoHs
<b>Supplied accessories</b>	User Guide Quick-Installation Sheet Video cables (1 x DVI Single Link, 1x DisplayPort) Mains cables (UK, European (CEBEC/KEMA), USA (UL/CSA; adaptor plug NEMA 5-15P), Chinese (CCC)) USB 2.0 cable External power supply
<b>Optional accessories</b>	NA
<b>QA software</b>	MediCal QAWeb
<b>Unit per pallet</b>	NA
<b>Pallet dimensions</b>	NA
<b>Warranty</b>	5 years
<b>Operating temperature</b>	0°C to 40°C (15°C to 32°C within specs)
<b>Storage temperature</b>	-20°C to 60°C
<b>Operating humidity</b>	8% - 70% (non-condensing)
<b>Storage humidity</b>	5% - 95% (non-condensing)
<b>Operating altitude</b>	3600 m
<b>Storage altitude</b>	7500 m
<b>Regulation caution</b>	 <p>Medical equipment with respect to electric shock, fire and mechanical hazards only in accordance with UL 60601-1. and CAN/CSA C22.2 NO. 601.1.</p>

\* In dark reading room conditions (0 lux)

# Warranty statement

## BARCO MID GENERAL WARRANTY TERMS AND CONDITIONS

Applicable from Jan 1st 2008

### ARTICLE 1: PRODUCT WARRANTY

Barco nv, Medical Imaging Division warrants that the equipment will be free of defects in workmanship or material for the warranty period or the specific period of a warranty extension program.

Hardware: Barco warrants that upon delivery hereunder the Products shall (i) conform to its specifications in effect at the date of delivery and (ii) be free from defects in material and workmanship (the "Warranties").

Software: Barco warrants that software written by Barco shall perform substantially in accordance with the specifications in effect at the date of delivery. Software is inherently susceptible to bugs and errors. Barco makes no warranties with respect to the software which is provided to Customer on an "as-is" basis and does not warrant uninterrupted or error-free operation of the Products.

Unless otherwise indicated in Barco's Product manual or in the agreement between Barco and Customer, the Warranty Period shall be

(a) Hardware: 12 months commencing on the Barco date of invoice.

(b) Software: 3 months commencing on the Barco date of invoice.

Notwithstanding the provisions of clause 2, repair and replacement of defects in material and/or workmanship under this warranty shall be accomplished in our works according to the terms and conditions as set forth hereafter:

1.1 Any claim under the Warranties must be notified to Barco in writing within 8 days from the date the defect or failure has been discovered or noticed the first time. The Customer, upon the occurrence of any equipment failure, shall contact Barco nv, MID customer support centre (or an authorised service centre) by telephone, fax or e-mail and shall provide the applicable customer support person with a complete description of the problem being encountered, including the model and serial number of the equipment in which the problem has arisen.

1.2 The customer support person shall diagnose the problem experienced by the Customer and shall advise the Customer on how to proceed. Customer support may ask to return the faulty equipment or faulty subassemblies to the Barco nv, or a MID customer support centre (or an authorised service centre) for repair activities. In no event shall Customer return a defective Product or part thereof to Barco without Barco's prior written approval.

In the event the return to Barco of the defective Product is authorized by Barco, Barco shall issue to Customer an RMA (Return Material Authorization).

The Customer shall apply for an RMA number to the closest Barco nv, MID Customer support centre (or an authorised service centre) as listed at [www.barcomedical.com](http://www.barcomedical.com), unless otherwise indicated.

The one-way cost of packing, transport and insurance related to shipping the alleged defective Product or part to Barco for repair or replacement shall be borne by Customer. The one-way cost of packing, transport and insurance related to shipping of the repaired or replacement Product or part to Customer shall be borne by Barco.

1.3 The Customer shall return, freight prepaid, the defective equipment or subassemblies in its original packaging with the assigned RMA number for repair to the Barco nv, MID Customer support centre (or an authorised service centre).

1.4 Replacement parts used shall be new or equivalent to new parts for the revision level of the equipment. A replacement LCD panel will be new or similar run time. The warranty period for the replacement parts will expire at the same moment as the original warranty period of the equipment. All parts replaced hereunder and returned to Barco nv, MID (or an authorised service centre) shall become the property of Barco nv, MID (or the authorised service centre).

1.5 The repaired equipment shall be returned to the Customer, by regular freight, at Barco nv, MID's charge.

1.6 Barco nv, MID will replace a product that fails within 90 days after shipment from Barco with a new one, ensuring color matching in the event of a diagnostic multi head display configuration.

Barco aims to ship the new replacement product within 2 business days after receipt of the Customers notification in writing.

### ARTICLE 2: CONDITIONS PRECEDENT FOR WARRANTY AND ITEMS EXCLUDED FROM WARRANTY

- A. The Warranties shall apply only to the extent the Products or any parts thereof have
- (i) been transported and stored at all times in the original packaging in the conditions as specified by Barco (such as covered and secure location, minimum and maximum temperature, maximum humidity, ...) or, in absence thereof, at least in conditions consistent with generally accepted practice for this type of products;
  - (ii) been handled at all times in accordance with Barco's instructions or, in absence thereof, at least with the care and caution consistent with generally accepted practice for this type of products;
  - (iii) been installed strictly in accordance with the instructions and directions given by Barco (if and to the extent the

- Products have not been installed by Barco or its authorized subcontractors);
- (iv) not been subject to any unauthorized access, alteration, modification or repair or attempts thereto;
  - (v) been at all times "normally used" for the intended purpose and operated in strict accordance with the operating instructions set forth in the operating manual of the Product in question and shall not have been otherwise misused, abused, damaged. For the purpose hereof, "normally used" shall mean a regular, ordinary and routine usage of the Product in question as intended and/or recommended by Barco;
  - (vi) been maintained at all times in accordance with Barco's instructions or, in absence thereof, at least with intervals and in a manner consistent with generally accepted practice for this type of products;
- B. The warranty described herein shall not include the following:
- 2.1 Any hardware or software item procured from a source other than Barco nv, MID or their official agent or distributor and integrated by Customer or a third party into Barco nv, MID supplied equipment.
  - 2.2 Any host configuration not explicitly supported by Barco nv, MID.
  - 2.3 All software installed on the system, whether they are acquired from Barco nv, MID or third party. An exception is made for software delivered by Barco nv, MID that would prove to be a cause for the mal functioning of the hardware covered under this Agreement.
  - 2.4 Normal wear and tear, use under circumstances exceeding specifications, such as use in dusty environment or under excessive temperature conditions, abuse, unauthorised repair or alteration, lack of proper configuration or maintenance, damaged or modified or removed serial number, cosmetic refurbishment.
  - 2.5 Repair or replacement of consumables<sup>1</sup> or specific parts that by definition are subject to wear and tear, including but not limited to:
    - a. CRT's, LCD panels
    - b. Backlights in diagnostic LCD displays, when the backlight run time<sup>2</sup> is beyond the Guaranteed Backlight Lifetime<sup>3</sup> of that model, when used at the Factory Calibrated Luminance<sup>4</sup>
      - Eg1.
        - A display is used at 8 hours/day; ie. +/-2.920hours/year
        - The display system is covered with a warranty period of 5 year
        - The Guaranteed Backlight Lifetime of that model is 17.000hours.
        - The Factory Calibrated Luminance cannot be achieved anymore after 4 year, thus corresponding to +/- 11.680hours
        - Result: the backlight replacement is performed under warranty
      - Eg2.
        - A display is used at 24 hours/day; ie. 8.760hours/year
        - The display system is covered with a warranty program of 5 year
        - The Guaranteed Backlight Lifetime of that model is 17.000hours.
        - The Factory Calibrated Luminance cannot be achieved anymore after 2.5years, thus corresponding to +/- 21.900hours
        - Result: the backlight replacement does not make part of the warranty coverage.
    - (1): 'Consumable' is a part that can be replaced by the user
    - (2): 'Backlight run time' is the total time that an image (including use of a screen saver) has been applied to the screen; this value can be consulted via the OSD buttons (On Screen Display).
    - (3): 'Guaranteed Backlight Lifetime' is the number of backlight hours during which a predefined luminance value, ie. the Factory Calibrated Luminance, will be provided by a specific diagnostic display model. This figure is printed on the Warranty certificate of each display or can be requested at your local Barco office.
    - (4): 'Factory Calibrated Luminance': is the typical luminance value that a specific diagnostic display model is calibrated at during the production process. This figure is printed on the Warranty certificate of each display or can be requested at your local Barco office.
  - c. Lamps, optical components in projectors
  - d. Replacement because of:
    - i. image retention as a result of:
      - not correctly using screen saver and/or Display Power Management System (DPMS) as explained in the user manual
      - prolonged operation of the display with a static image on the same screen area
    - ii. a number of missing pixels that is lower than the total allowable number as mentioned in the product specifications.



- iii. difference in color temperature that is lower than the total allowable difference as mentioned in the product specifications
- iv. difference in color temperature as a result of not using all displays of a multi-head configuration at the same rate.
- v. LCD luminance uniformity that is in within the product specifications or luminance uniformity performance that is inherent to LCD technology.

Barco nv, MID does not warrant a minimum life time nor a performance of any of the consumables.

2.6 Replacement of moveable parts such as power cords, remote controls, ...

2.7 Any product disassembly and installation costs at the operation site, travel expenses and travel time to and from the operation site for the personnel in charge of the repair works and transport charges.

2.8 Any failures resulting from an accident, negligence (such as but not limited to removing or deleting system files & licensed software product files), misuse, circuit failure or any change, damage due to fire, water, thunder or lightning, power failure or fluctuation, disruption of communication lines or due to force majeure, or any reason foreign to the equipment.

2.9 Any specific services or procedures, asked for by Customer, related to verification of repaired equipment.

2.10 The evaluation cost in case of a returned product deemed functional is not covered under warranty and will be charged to the Customer, at Barco's sole discretion, based on a case-by-case evaluation.

### **ARTICLE 3: OBLIGATIONS OF THE CUSTOMER**

Customer hereby assumes the following obligations as partial consideration for Barco nv, MID performance of its requirements under the warranty condition; failure by Customer to meet its obligations under this paragraph shall excuse Barco nv, MID's performance hereunder:

3.1 Customer shall not expose Barco nv, MID personnel to any unsafe working conditions.

3.2 Repairs to equipment under warranty resulting from improper maintenance or repair performed by the Customer, or its officers, agents, employees, or representatives, shall be borne by the Customer at its additional cost and expense.

3.3 The Customer is responsible for installing the Barco nv, MID equipment in an environment for which it was intended. If there is an indication that the equipment was used – even temporary – outside its specifications, Barco nv, MID is entitled not to perform warranty repairs and terminate the warranty agreement. Any actions that have been taken by Barco nv, MID in this respect, may be invoiced to the Customer at normal pricing.

### **ARTICLE 4: REMEDIES UNDER WARRANTY**

(a) Hardware: If during the Warranty Period a Product or any part thereof, fails to meet any of the Warranties then, upon Customer's request, Barco shall, at its sole option and cost, promptly and within 20 working days, either: (i) repair or correct the Product or part in question; or (ii) replace the Product or supply part(s) or component(s) according to the terms and conditions contained in article 1. A replacement part shall be at least functionally equivalent to the original part. The replaced Product, parts and/or components shall become the property of Barco and shall, at Barco's request, be returned by Customer to Barco at Barco's cost.

(b) Software: Barco's sole obligation shall be to rectify substantial malfunctions of the software (to the extent technically reasonably possible) by amending the software or supplying an alternative version of the software.

(c) The repair or replacement under the Warranties covers the cost of material and labor.

### **ARTICLE 5: DISCLAIMER OF WARRANTIES**

BARCO NV, MID DISCLAIMS ALL WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING ALL IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

### **ARTICLE 6: LIMITATION OF LIABILITY**

The remedies specified in these Terms and Conditions shall constitute Customer's sole and exclusive remedy and Barco's sole and exclusive liability for Barco's breach of the Warranties hereunder.

BARCO NV, MID SHALL NOT UNDER ANY CIRCUMSTANCES BE LIABLE TO CUSTOMER OR ANY THIRD PARTY FOR DIRECT, INDIRECT, INCIDENTAL, SPECIAL OR CONSEQUENTIAL DAMAGES, SUCH AS BUT NOT LIMITED TO, DAMAGE TO OR LOSS OF TANGIBLE OR INTANGIBLE PROPERTY OR EQUIPMENT, LOSS OF PROFITS OR REVENUES, COST OF CAPITAL, COST OF PURCHASE OF REPLACEMENT GOODS, OR CLAIMS OF USER FOR SERVICE INTERRUPTIONS. THE LIABILITY OF BARCO NV, MID FOR MANUFACTURING, SALE, DELIVERY, RESALE, INSTALLATION, OPERATION OR SUITABILITY FOR USE OF ANY PRODUCTS OR SERVICES COVERED BY OR FURNISHED UNDER THIS WARRANTY CONDITION, WHETHER ARISING OUT OF CONTRACT, NEGLIGENCE, STRICT TORT, WARRANTY OR OTHERWISE, SHALL NOT EXCEED THE PRICE OF THE ITEM OR ITEMS OF GOODS OR SERVICES UPON WHICH SUCH LIABILITY IS BASED.

### **ARTICLE 7: FORCE MAJEURE**

Either party shall be released from performance of its obligations under this agreement to the extent, and for so long as, the performance of this agreement is impeded by reason of force majeure. For the purposes of this clause the

expression "force majeure" means, but shall not be limited to, industrial dispute, fire, mobilisation, requisition, embargo, currency transfer prohibitions, insurrection, lack of means of transport, restrictions of the use of energy, and generally any circumstances which are beyond the control of the parties and hinder performance by one party of his obligations.

#### **ARTICLE 8: GENERAL**

8.1 Customer acknowledges its understanding that all software and electronic devices, including Barco nv, MID products are subject to possible error, mechanical or electrical failure, and should not be relied upon in inappropriate applications or without proper backup and/or other safety precautions whenever personal injury or property damage may result from failure or error of the product.

8.2 Barco nv, MID shall not be responsible for machine failure and/or its failure to render service or maintenance due to causes beyond its reasonable control.



**B410582-08**

March 2012



**Barco n.v.**  
President Kennedypark 35  
8500 Kortrijk  
Belgium  
[www.barco.com](http://www.barco.com)